


**The ez
Guide**

Model # UVC 3069

How to setup:

99.9% sterilizing rate

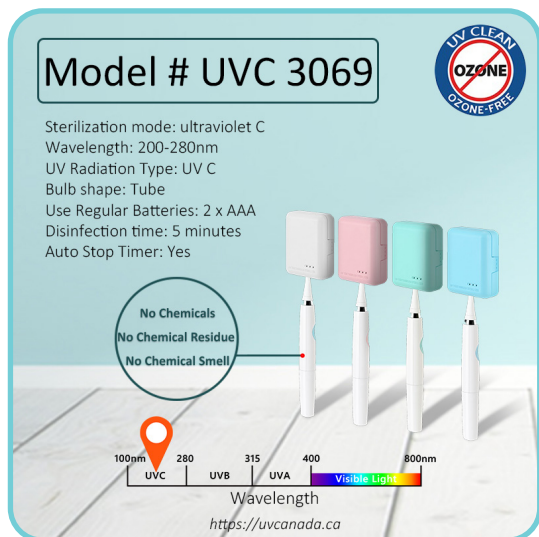

Step 1: Choose the appropriate fixed soft glue in the packaging that fits the width of the neck of the toothbrush

Step 2: Open the battery cover located at the bottom of the unit and install 2x “AAA” alkaline battery

Step 3: Open the cover of the toothbrush sterilizer

Step 4: Place the head of the toothbrush side up

Step 5: Close the cover for 5 minutes for eradication

Note: The unit will stop working if the cover was not closed properly

Benefits

Suggestions

- UV-C Fast acting and efficient sterilization
- Unique and simple mode of operation
- 2 x AAA batteries
- Eliminate a variety of Pathogens
- Free from chemical agents
- User-friendly
- Ozone-Free

- Best Use: Toothbrush

Caution and Warnings

- Do not look at the light source with naked eye; avoid light from falling on skin
- Avoid exposure to UV rays, it can be harmful to the eyes and skin
- Do not use this Sterilizer for other than its intended purpose, as described in suggestions
- Equipment should be protected from the elements and from temperatures below freezing
- Do not operate the UV Sterilizer if any damage to the unit is visible or suspected
- Do not alter the device
- This equipment is intended for indoor use only

Maintenance

Keep the UV Sterilizer clean you can gentle wipe it with soft cloth